



Brushlink[®]

YOUR PERSONAL DENTAL COACH

INSERTING BATTERIES

1

Remove the lower blue-coloured casing.



2

Insert the two LR44 batteries provided. Ensure that the batteries have been inserted as per the image below:



3

Place the white-coloured upper casing on top of the lower casing to close the device.




ATTACHING BRUSHLINK®



Align Brushlink on the back of the toothbrush handle, with the light facing up.

BAND SIZES

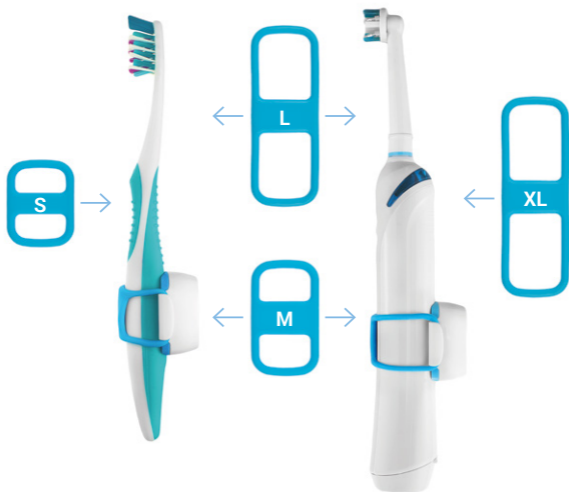

Wrap the **silicon band** around the toothbrush handle to secure Brushlink into place.

Manual toothbrushes use the small or medium size blue band. **Electric toothbrushes** use the medium, large or X-large.

DOWNLOADING THE BRUSHLINK® APP

The Brushlink app is available on both Android and Apple iOS mobile devices.

Click on the icons below to download the app for your smartphone:



SWITCHING BRUSHLINK® ON



Press the **blue button** to activate your device and connect it to the Brushlink app.

ENSURE THE:

- **App is open** on your smartphone and
- **Bluetooth™ is turned on** to allow the app to automatically connect to your device.

LIGHTS ON THE BRUSHLINK®



GREEN LIGHT

The Brushlink device is switched on.



BLUE LIGHT


Brushlink is connected to the app, via Bluetooth.




RED LIGHT

Deleting data from the Brushlink device or low battery.


BRUSHLINK® APP SET UP

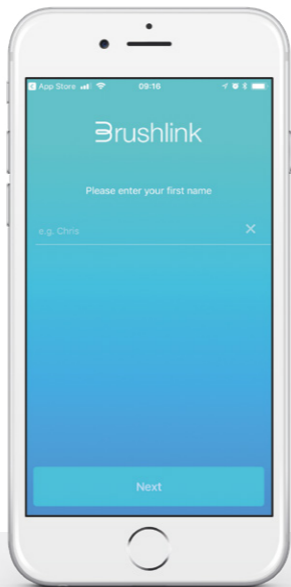


Fill in your personal details.

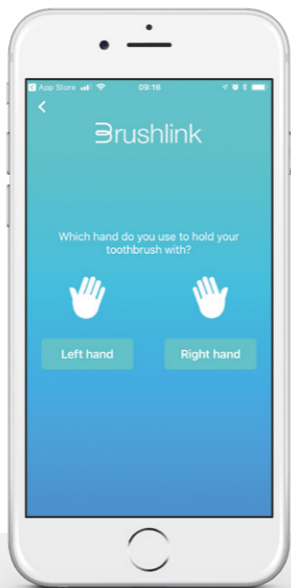


Choose left or right-handed.

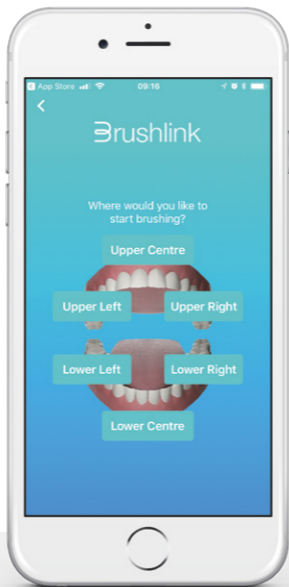

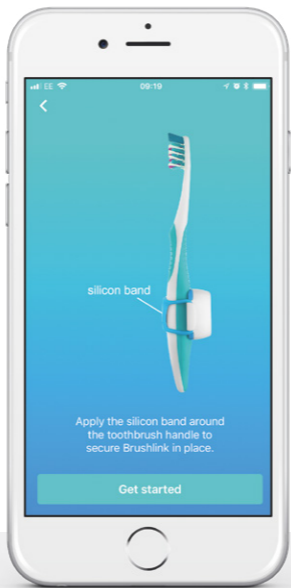


Choose an area to start brushing.

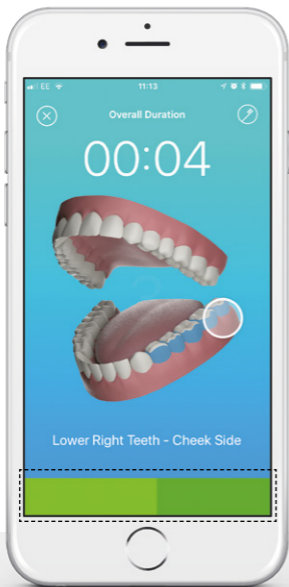


Attach the Brushlink to your toothbrush (if you have not done so).




The bar at the bottom of the screen shows how **well** you are brushing.

-  Correct angle
-  Improvement needed
-  Incorrect angle

Brush your teeth following the app.

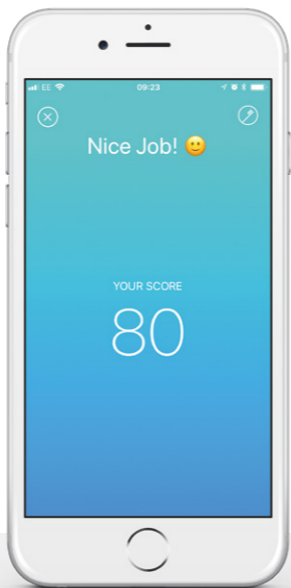

Once you are finished you will be marked out of 100.

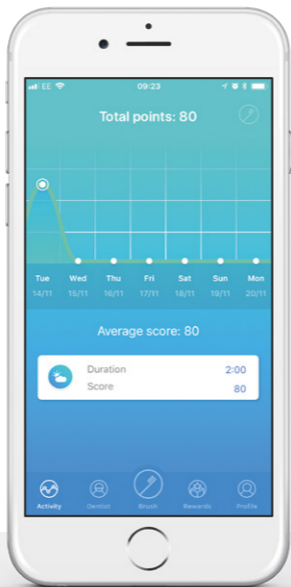

Click on Activity to view your score.

HOW TO BRUSH EFFECTIVELY



02:00

Make sure you brush for at
least two minutes.

To get the best brushing angle,
the bar should **always be green.**

HOW TO BRUSH EFFECTIVELY



Ensure your toothbrush is at a 45° angle to your teeth and gums.



The background is a solid blue color with various abstract light effects. On the left, there are curved, glowing lines that resemble a stylized 'B' or a signal waveform. Scattered throughout the background are numerous small, bright white and light blue dots, some of which are larger and more prominent, creating a bokeh effect. The overall aesthetic is clean, modern, and tech-oriented.

brushlink.com

© DNC Technologies Limited. All rights reserved.